

Inspection Report

Designation	: Standard droplet Sample (Calibration silhouette)
Specification	: Standard View
Application	: Calibration for contact angle meter
Material	: Chromium sputtering prints on Soda-lime glass
Manufacturer	: Kyowa Interface Science Co., Ltd.
Serial No.	: DSD9999
Date of Inspection	: February 6, 2024
Environment	: temperature: 23.1℃ humidity: 51%

Acceptable standard There must not be obstacle for the usage.

The diameter of the globe is measured three times both in horizontal and vertical directions, and variations of all the measured values must be smaller than 0.003mm and that all the measured values must be within ± 0.005 mm of the indicated values.

Indicated value 2.999 mm

	1	2	3
Diameter in horizontal(mm)	3.000	3.000	3.000
Diameter in vertical(mm)	2.998	2.998	2.998

Difference between the largest and smallest measured value 0.002 mm

Largest difference between the indicated value and measured values 0.001 mm

Conclusion Acceptable

Inspector and standards for inspection.

Inspector	: Pronics Co.,Ltd
Item	: CNC Video Measuring System NEXIV
Model	: VMR-6555
Manufacturer	: NIKON SOLUTIONS CO., LTD.
Serial No.	: 1322533

Approved by :



Kyowa Interface Science Co., Ltd.

Calibration Certificate

Designation	: Standard droplet sample (Calibration silhouette)
Specification	: Standard View
Manufacturer	: Kyowa Interface Science Co., Ltd.
Serial No.	: DSD9999
Date of Inspection	: February 6, 2024
Environment	: Temperature : 23.1℃ Humidity : 51%


This is to certify that the above instrument has been calibrated and inspected on the basis of the standard of Kyowa Interface Science Co., Ltd., and the results are conformity with the Inspection Reports enclosed herein.

The standard tool used for calibration and inspection has been calibrated on the system of the calibration and inspection of Kyowa Interface Science Co., Ltd., and has been traceable to the national standard of Japan as the Traceability System enclosed herein.

Inspector and standard for inspection.

Inspector	: Pronics Co.,Ltd
Item	: CNC Video Measuring System NEXIV
Model	: VMR-6555
Manufacturer	: NIKON SOLUTIONS CO., LTD.
Serial No.	: 1322533




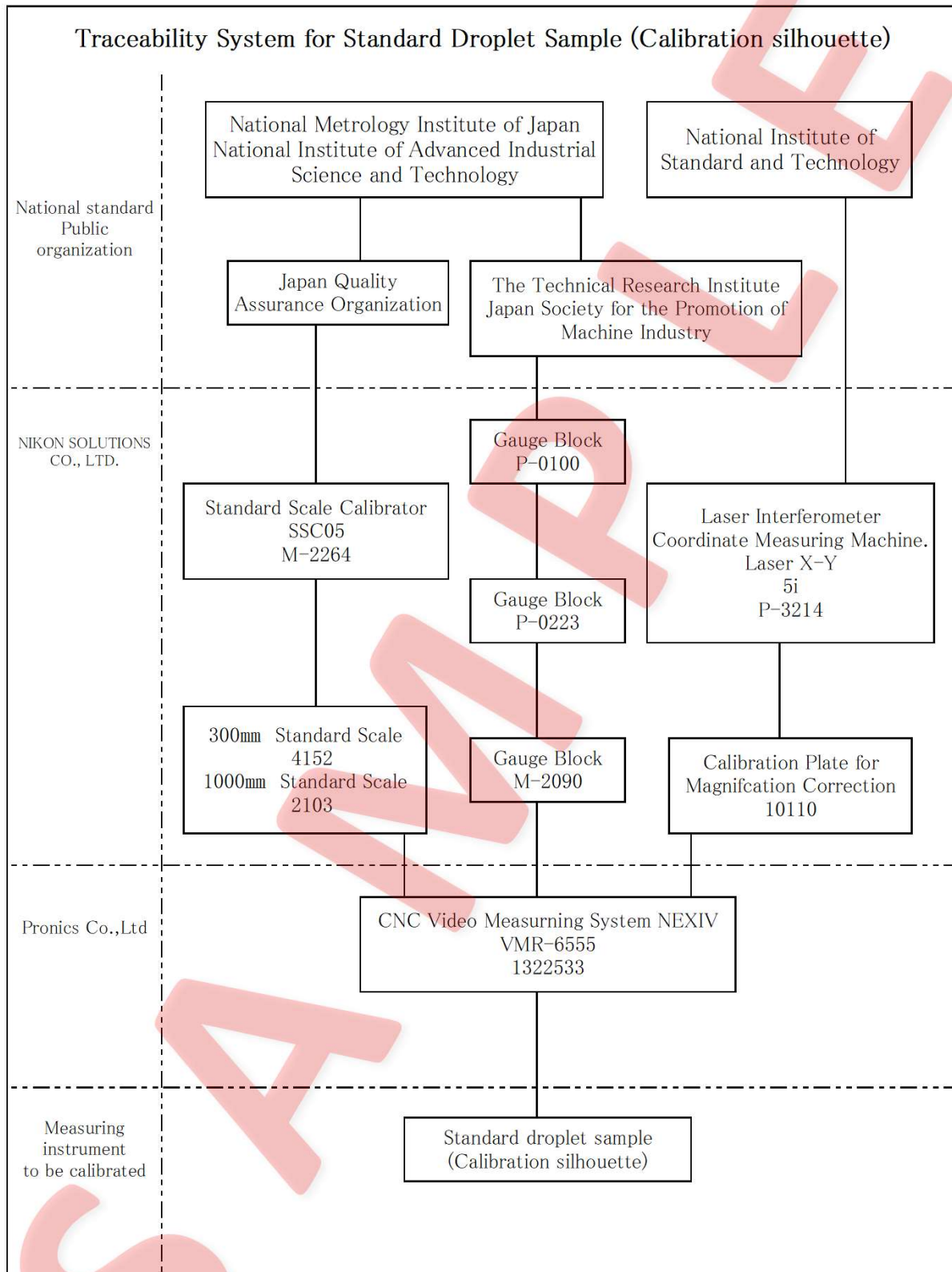

Inspection Report

Designation	: Standard droplet sample (Droplet silhouettes)
Specification	: Standard View
Application	: Calibration for contact angle meter
Material	: Chromium sputtering prints on Soda-lime glass
Manufacturer	: Kyowa Interface Science Co., Ltd.
Serial No.	: DSD9999
Date of Inspection	: February 6, 2024
Inspection environment	: Temperature : 23.1℃ Humidity : 51%

Acceptable standard: There must not be obstacle for the usage.
The major axis (2r) and the height (h) of the droplet silhouette are measured. Then, contact angle is calculated from those measurement data and it must be within +/- 0.5 degree against the nominal value.

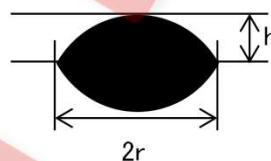

Nominal value	Major axis (2r)	Height (h)	Contact angle (θ)	Conclusion
5.0	4.923	0.107	4.98	Acceptable
60.0	3.165	0.913	59.96	Acceptable
108.0	2.087	1.434	107.91	Acceptable

Inspector and standard for inspection.

Inspector	: Pronics Co.,Ltd
Item	: CNC Video Measuring System NEXIV
Model	: VMR-6555
Manufacturer	: NIKON SOLUTIONS CO., LTD.
Serial No.	: 1322533

Approved by :



Kyowa Interface Science Co., Ltd.

Calibration Certificate

Item	: Standard droplet sample (Droplet silhouettes)
Specification	: Standard View
Manufacturer	: Kyowa Interface Science Co., Ltd.
Serial No.	: DSD9999
Date of Inspection	: February 6, 2024
Environment	: Temperature : 23.1°C Humidity : 51%

This is to certify that the above instrument has been calibrated and inspected on the basis of the standard of Kyowa Interface Science Co., Ltd., and the results are conformity with the Inspection Reports enclosed herein.

The standard tool used for calibration and inspection has been calibrated on the system of the calibration and inspection of Kyowa Interface Science Co., Ltd., and has been traceable to the national standard of Japan as the Traceability System enclosed herein.

Inspector and standard for inspection.

Inspector	: Pronics Co.,Ltd
Item	: CNC Video Measuring System NEXIV
Model	: VMR-6555
Manufacturer	: NIKON SOLUTIONS CO., LTD.
Serial No.	: 1322533



